

DATASHEET

FFP MASKS

PROTECTION AGAINST DUST, MIST & FUMES



SMART SERIES – SMART SOLO MASKS

FFP2 NR D



2495 with Ventex®-valve

FFP3 NR D




2595 IonicAir® with Ventex®-valve

CHARACTERISTICS



ActivForm®

Automatically fits to the face.
No manual adjustments by the user are necessary.

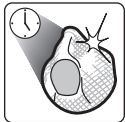


DuraMesh®

Masks have a strong and durable structure.




Ventex®-valve

Starts to open even at low exhalation pressure and significantly reduces moisture and heat inside the mask.

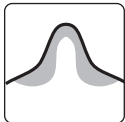

Nose seal

The flexible nose seal improves fit and provides optimum wearing comfort.

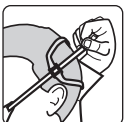


1-Strap

Designed for quick and simple fitting and removal of the mask, even when wearing gloves.



Flexi-Wings

Evenly distributes the strap-force to ensure a safe fit.

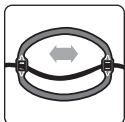

Head Harness

The adjustable head harness ensures correct positioning and optimal wearer comfort.



Dolomite clogging test

Masks have passed the Dolomite clogging test.
Better breathing resistance for longer.



100% PVC-FREE

All Moldex products and packaging are completely free from PVC.

CERTIFICATION

The Moldex Smart Solo FFP-masks meet the requirements of EN 149:2001+ A1:2009. The products are CE-marked in accordance with the requirements of EU regulation (EU)2016/425. The IFA (0121) in St. Augustin (Germany) is responsible for type examination (Module B) and monitoring of production (Module D). The products are manufactured in an ISO 9001 certified plant.

MATERIALS

Filter Layer: Polypropylene

Inner Shell, DuraMesh®: Polypropylene, Ethylene vinyl acetate (EVA)

Nose Seal, Head Harness: Polyethylene

Ventex®-valve: Natural Rubber

Head Strap: Polyester, Lycra

WEIGHT

2495: 25 g

2595 Ionic Air®: 26 g

AREAS OF USE

Level	WEL	Hazard type
FFP2	10 x	Hazardous fine dusts, water and oil based mists/aerosols, Biological agents of risk group 2
FFP3	20 x	Harmful and carcinogenic dusts, water and oil based mists/aerosols, biological agents of risk group 2 and 3, CMR-substances

(WEL = Workplace Exposure Limit)

NR (non reusable) = Single use. Comfortable and durable throughout the whole shift

DATASHEET

FFP MASKS

PROTECTION AGAINST DUST, MIST & FUMES



TESTING ACCORDING TO EN 149:2001 + A1:2009

Total inward leakage

Ten test subjects perform a variety of exercises. During the exercises the amount of test aerosol that penetrates the filter, face seal and valve are sampled. The total inward leakage of 8 out of 10 test subjects shall not exceed the following levels:

Category	max. total inward leakage
FFP2	8 %
FFP3	2 %

The filter penetration after loading the filter with 120 mg paraffin oil according to EN 149:2001 + A1:2009 shall not exceed the following levels:

Category	max. Filter penetration
FFP2	6 %
FFP3	1 %

Flammability

4 respirators are passed through a 800°C (+/- 50°C) flame with a speed of 6 cm/s. After passing through the flame the respirator has to self-extinguish.

Breathing Resistance


The breathing resistance produced by the filter of the respirator is tested at an airflow of 30 l/min and 95 l/min.

Category	max. breathing resistance according to EN 149	
	30 l / min	95 l / min
FFP2	0,7 mbar	2,4 mbar
FFP3	1,0 mbar	3,0 mbar

INSTRUCTIONS FOR USE

- The user has to be trained and instructed in wearing the mask.
- FFP masks do not protect against gases and vapours.
- The oxygen concentration of the ambient atmosphere should be at least 19,5 % Volume.
- These respirators may not be used if the concentration type, and properties of contaminants in the ambient atmosphere are unknown or at dangerous levels.
- Respirators should be disposed if damaged, if the breathing resistance becomes high due to clogging, or at the end of a shift.
- Never tamper with, alter or modify the respirator.


INSTRUCTIONS FOR FITTING

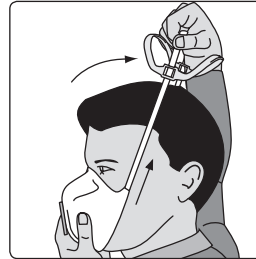

1. Place respirator on chin and pull head strap on the head harness over the crown of the head.



2. Place head harness on the back of head.

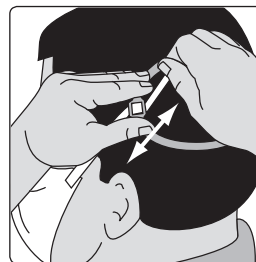

3. Ensure the respirator fits securely and comfortably.
If necessary, adjust strap by pulling it at either side of the head harness.

INFO

For help on selection and training please contact us. We offer a wide range of training packages and support material.

MOLDEX/METRIC AG & Co. KG
Unit 9, Glaisdale Point
Off Glaisdale Drive, Bilborough
Nottingham. NG8 4GP

Tel.: +44 (0)115 985 4288
Fax: +44 (0)115 985 4211
info@uk.moldex-europe.com
www.moldex-europe.com